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Dr. Mayer
Sensitive
Pressure
GTS2099

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Dr. Mayer Sensitive Pressure
GTS2099

BetterClean • Better Dental Care

Before use, please carefully read the instructions. Retain the instructions for future reference.

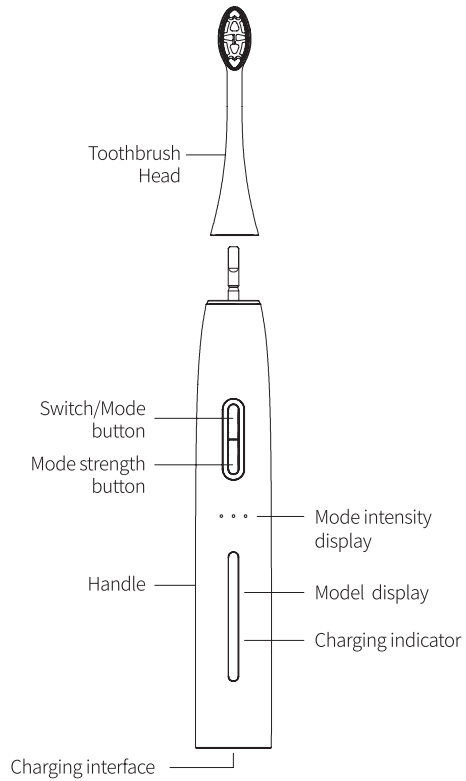
DEAR CUSTOMER,

Thank you for choosing Dr. Mayer products. Please read this instruction before use and use it in a proper way so as to achieve ideal effects.

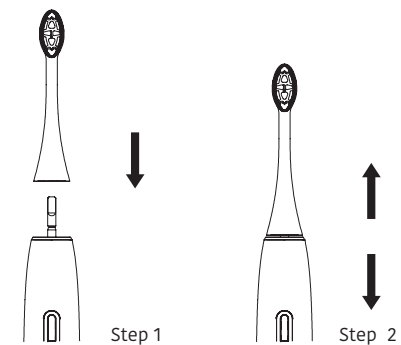
IMPORTANT DAFEGUARDS

- There are some people may not adapt (such as stimulation and feeling numb) when using for the first time, please don't worry, this is the normal phenomenon, this phenomenon will fade out after several times using, and you will be adapt to this toothbrush in one week.
- Do not open, repair or alter the product without authorization it.
- Do not knock or drop the product, if product is defective, please stop using.
- Please check if the local voltage is suitable with the marked voltage before connecting the power supply. This product can only use the special toothbrush head. If your oral or gum had surgery during the past two months, please consult your dentist before use.
- If your teeth is bleeding after using the toothbrush or bleeding continues more than one week, please consult your dentist.
- This product is designed for teeth, gums and tongue cleaning, not for other purpose.
- Children should be supervised to ensure proper usage and make sure they do not play with the toothbrush. Do not place in high temperature condition or directly shining under the sunlight.
- Please keep the charging port dry while charging and use the original charging cable!

COMPONENT NAME

INSTRUCTION

PRODUCT INFORMATION

Product name: Sonic electric toothbrush
Model: Dr. Mayer Sensitive Pressure

Input voltage: 5V
Input power: 2W
Product material: ABS / TPE / Nylon
Toothbrush size: 30 × 28 × 250 mm
Charging time: ≤ 8 hours
Working temperature: 0-50 °C
Orgin: China

FEATURES

1. Clean oral effectively with 42,000 times per minute.
2. The product has 5 modes: white, clean, sensitive, polish, massage. Each brushing mode has three intensity levels to fit different needs.
3. Pressure sensor feedback: if too much pressure is applied, the vibration will be reduced and all indicator lights will flash.
4. Safe and durable lithium battery runs for about 60 days on a full charge (in the strongest mode protection), wireless charging, with a power short-age reminder.
5. Anti-splash function: when the product is turned on, the vibration will be low for 3 seconds, preventing toothpaste from splashing. Press the button within 3 seconds to exit the antisplash function.
6. The product is anti-skid, easy to hold and ergonomically designed.
7. Waterproof class at IPX7.
8. Smart mode memory, 2 minutes auto stop and 30 seconds timer reminder for moving help you get good brushing habits.

USING THE TOOTHBRUSH

1. Attach the brush head and make sure it is locked in place. Wet the bristles and apply appropriate amount of toothpaste.
2. Gently put the brush head on your teeth.
3. Gently move the brush head across each section of the teeth for 30 seconds, including front, back and chewing surfaces.
4. The multifunction button needs continued operation in 3 seconds for first mode switching, press the strength key to select strength, otherwise, it will automatically stop working.
5. Switch from 5 modes (white, clean, sensitive, polish, massage) and 3 intensity levels (low, medium, high) according to your own preference.
6. Move to next section when the auto reminder pause briefly every 30 seconds. After finishing all 4 sections, it will automatically stop after 2 minutes.

CHARGING INSTRUCTION

1. If only one indicator light is on after use, or if the product cannot be turned on normally, it indicates that the product needs to be recharged. Please use a compatible charging cable to charge the product.
2. When charging, the mode indicator lights repeatedly flash from bottom to top. When charging is complete, all indicator lights are off.
3. For the first time usage or when battery is undercharged for long period of time, it is normal the indicator light does not light for initial 5 minutes.

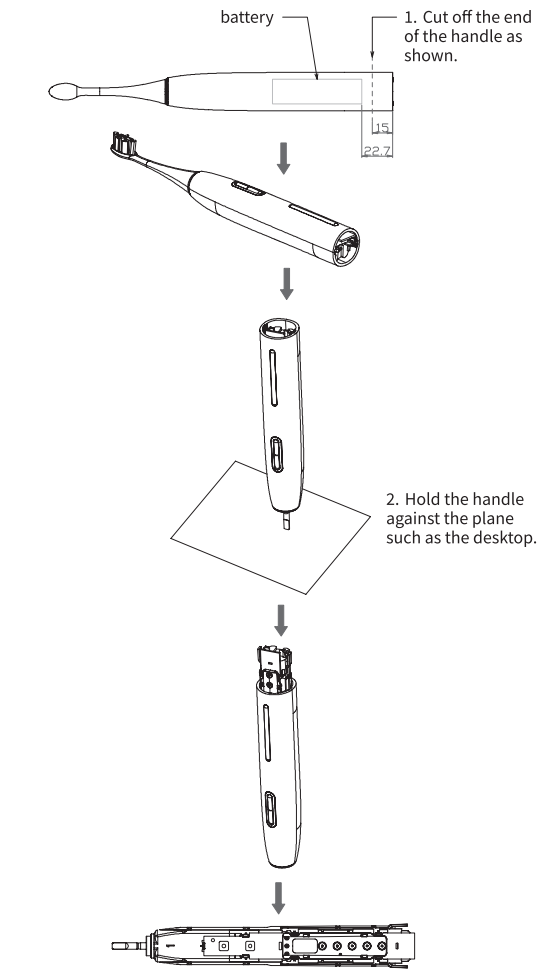
ENVIRONMENTAL PROTECTION RECYCLING

METHOD OF RECHARGEABLE BATTERY

This product contains lithium battery. Before being discarded, please take out the battery and recycle or dispose according to the local waste a compatible charging cable to charge the disposal management regulations. If in doubt, consult your local waste management department.

When removing the battery, you need the following tools: pliers, scissors.

Notices: when you decide to discard the toothbrush, you need to take out the battery. Make sure the battery is discharged before removing. When using the tools, please observe the basic safety rules and take protective measures to avoid injury.

3. Remove the battery from the holder after taking out the bracket.

TROUBLE SHOOTING

Case	Possible cause	Solution
The work time is short.	1. The handle is not in the right place which causing the charging failure.	1. Ensure to put it in a right way when charging.
	2. The battery maybe charge fail.	2. To extend the charging time.
The product doesn't work.	1. Checking if the product is full charged.	1. Please charge the battery.
		2. Please Contact the seago customer service to deal with it.